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emmi[®]-dent Metallic

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Emmi Ultrasonic GmbH • Gerauer Strasse 34
64546 Moerfelden-Walldorf • Germany
Tel: +49 (0)6105 - 40 67 94
Fax: 49 (0)6105 - 40 67 81

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SUPERIOR DENTAL HYGIENE

emmi®-dent - The original for innovative dental cleaning with 100% original ultrasound

We thank you for your purchase of an emmi®-dent Metallic Ultrasound device for your dental cleaning and oral hygiene („Ultrasound toothbrush“). Please read carefully this user manual for optimal results.

Content emmi®-dent Metallic Basic-Set

1 emmi®- dent Metallic hand piece,
1 emmi®- dent attachment for adults,
1 charging station,
1 tube ultrasonic toothpaste fresh(75ml)

Loading and service life per load

Important!

For each charging, connect the plug (type C) of the charging station into the socket (100 - 240 V), insert the hand-piece into the dock. When the charging station is connected properly, the green LED light on the charging station is ON. Charge the hand-piece once 24 hours for the first time, to activate the maximum performance of the battery. **The blue LED does not necessarily turn off after the 24 hours of charging.** If you use the device 6 minutes daily, you can use it with the battery for approx. 12 days. Charge it then as needed.

Before first use

Insert the ultrasound brush head (2-pin plug system) vertically in the hand-piece. Also remove the Ultrasound brush head in vertical direction.

Switching on and off

Switching on: Switch on the unit by pressing the pushbutton (for approx. 1 second). This emmi®-dent Metallic offers the function of switching the vibration motor and the position change signal (signal tone/beep tone) on and off so that you can decide whether you want to perform oral hygiene with or without the vibration motor and the signal tone. To switch off the vibration motor and the signal tone, please keep the On/Off button pressed for about 9 seconds. The LED lights up red four times. To switch the vibration motor and the signal tone back on, press and hold the On/Off button again for about 9 seconds. The LED lights up green four times.

Switch off/Automatic switch off: Switch off the unit by pressing the push-button (for approx. 1 second). The unit switches off automatically after approx. 6 minutes (Battery Save Function). The device is also provided with an electronic sensor. If a wrong or defective Ultrasound brush head is used, or no brush head has been inserted, the unit will also switch off automatically and the green LED goes off after 1 second. If you have forgotten the attachment on the handpiece for your trip, the device will also switch off after approx. 6 minutes to prevent discharging (Travel Safe function).

LED Display

1. ON:

You switch on the device by pressing the ON/OFF button. The green LED lights up (power indicator).

2. LED green - LIGHTS:

The green LED indicates that the emmi®-dent ultrasonic toothbrush works properly. Visible only in operation, not during charging.

3. LED blue - CHARGE:

The LED display lights up blue when the handset with the attached attachment is on the charging station and the charging station is properly supplied with 100 - 240V mains voltage. We use a type C plug.

4. LED red - CHECK/ DISRUPTION:

If the battery is empty or if there is a fault, the red LED display lights up. The device switches off, the red LED flashes for a few seconds and goes out.

Application of the emmi®-dent Metallic 2.0 with 100% original emmi®-dent ultrasound.

After switching on the emmi®-dent, you will feel the discreet vibration motor and a quiet position change signal every 6 seconds. The actual cleansing process is not affected by this. Ultrasound in the range of the multi-frequency bandwidth used by us is basically vibration-free and noiseless. Both functions serve to give the user feedback that the device is running. Also, a low position change signal tone is emitted every 6 seconds and the motor stops briefly to give you a haptic signal. (Beep can be switched off manually).

Daily use

For your dental hygiene, we recommend you to clean your teeth with the 100% Original Ultrasound of your emmi®-dent Metallic twice a day for approx. 3 min.

1. Moisten the bristles of the Ultrasound brush shortly under running water.
2. Apply a double pea-sized quantity of emmi®-dent Ultrasound Toothpaste on the bristles of your Ultrasound brush head. Only the emmi®-dent Ultrasound Toothpaste guarantees the Original 100% emmi®-dent Ultrasound function. For more details on the emmi®-dent Ultrasound Toothpaste, please read under the separate section „emmi®-dent Ultrasound Toothpaste“ in this user manual.
3. Switch on the device and bray the Ultrasound Toothpaste uniformly on the inside and outside surfaces as well as on the chewing sides of your teeth.
4. The bristles of the Ultrasound brush head embrace 2 – 3 teeth, according to the size of your teeth. Maintain the bristles gently and without pressure motionless on the inside, outside and on top (chewing sides) of the teeth for approx. 6 seconds each time. The actual duration depends on the amount of impurities on the teeth (plaque, nicotine, wine stains, etc.). In case of acute problems, such as gingivitis or parodontitis, leave the bristles for approx. 15 to 20 seconds on the affected sites. This can support the healing process. For more details please see our study at www.emmi-dent.com/info/studies. A longer application does not harm your teeth, because the emmi®-dent Metallic implements a gentle Ultrasound. When you have finished your dental hygiene, move the bristles various times over your teeth, in order to remove any further adherences or impurities. Then rinse your mouth with water or an oral rinse.
5. After cleansing, wash the bristles and hand-piece under running water. Switch off the device and dry the hand piece. Please note that the bristles are free and clean from foodstuff or other impurities. Now you can place the device with the integrated Ultrasound brush head back on the charging dock.

WITH 100% ULTRASOUND

1.

2.

3.

4.

5.

Gingivitis, periodontitis and aphthae in the mouth

In acute cases, it is recommended to hold the bristles gently for 30 to 60 seconds on the affected sides. This will result in a quicker and better healing. The Original emmi®-dent Ultrasound is harmless and can also be used for a longer period of time. We always recommend that you discuss your gum inflammations, periodontitis and aphthae with your dentist.

Bad Breathe

Bad breath (halitosis) is caused by living bacteria in the mouth. Since ultrasound basically reduces bacteria, bad breath can be prevented by using emmi®-dent.

emmi®-dent ultrasonic toothpaste


emmi®-dent ultrasonic toothpaste is available in five different versions: emmi®-dent Fresh, emmi®-dent Whitening, emmi®-dent Kids, emmi®-dent Nature and emmi®-dent Nature Calcium. Our Nature ultrasound toothpastes do not contain fluorides. All emmi®-dent ultrasonic toothpastes are gluten-free and contain no nanomaterials. To ensure the effectiveness of ultrasonic cleaning, microbubbles were added to all toothpastes from our company in an additional production step. For more information, please visit www.emmi-dent.com/info/studies.

SUPERIOR DENTAL HYGIENE

To transmit ultrasound, a transmission medium is required. It is a well-known concept in medical implementations. In this domain, creams and gels are used as adequate transmission medium. For your dental hygiene with Original emmi®-dent Ultrasound, you use a special toothpaste (emmi®-dent Ultrasound Toothpaste). This toothpaste stands out for its ultrasound feature, because emmi®-dent Ultrasound Toothpaste is the only medium in which millions of micro-bubbles are located. Conventional toothpastes in the market only generate a little amount of micro-bubbles. Always consider an emmi®-dent Ultrasound Toothpaste for your dental care with the emmi®-dent Ultrasound devices.

Replacement ultrasonic attachments

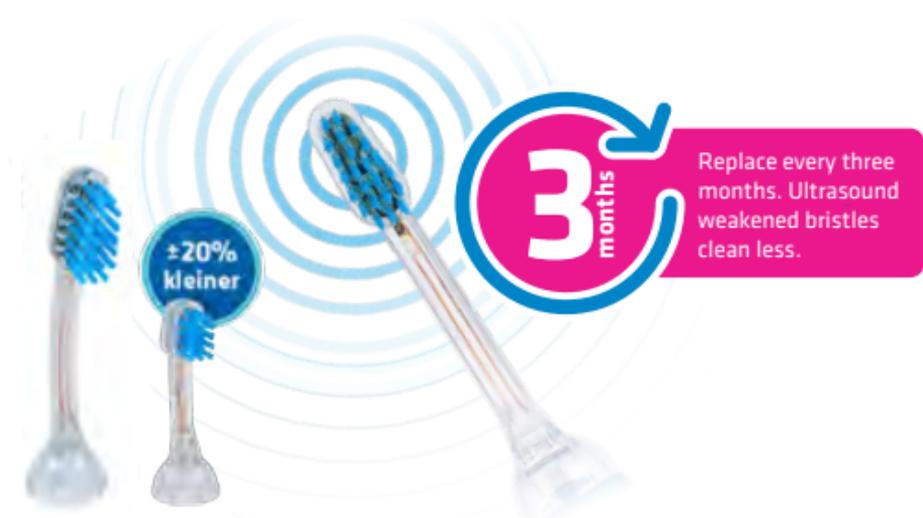The following ultrasonic attachments are available:

- emmi®-dent E2/E2 Wave (adults),
- emmi®-dent K2/K4 (children or women),
- emmi®-dent SB2 clasp carrier).

Further information is available at:

www.emmi-dent.com/equipment/replacement-brushes-for-metallic-professional/

Please take into account that according to dentists, the Ultrasound brush head should not be used by more than one person. One and the same Ultrasound brush head must always be used by one and the same person. Spare brush heads can be purchased at the locations where you have bought your emmi®-dent Metallic. The ultrasonic attachments should be changed after about three months, because the bristles are stressed by the ultrasound, which worsens the ultrasonic transmission and thus the result of the tooth cleaning.

How it works

Our 100% Original Ultrasound devices are the most innovative and effective method for a perfect dental care and oral hygiene (nano-cleansing). The emmi®-dent ultrasonic devices have been developed by German engineers according to the latest scientific findings and do not work like standard toothbrushes by means of mechanical friction and abrasive substances in the toothpaste (macro tooth cleaning), but use the properties of modern ultrasonic technology,

without pressure, without mechanical movements and thus without abrasion. It has been applied for brain, oftalmologic, pregnancy and stomach research. Ultrasound for dental and mouth hygiene is easy to use, whether we are young or old, because it is absolutely harmless and innocuous. Also carriers of pacemakers can use it without any risk. Research has proven that ultrasound of up to 50 watt is harmless for the human body. The emmi®-dent Metallic works with only 0,2 watt. The performance and safety of Ultrasound waves have been reviewed and confirmed by dentists, implantologists and studies. Many years of experience have never reported any negative effects, also for crowns, inlays, implants or veneers. The ultrasound waves are created inside an piezo-ceramic disk (ultrasonic transducer) which is patented. This chip is embedded inside the brush head and transmits the waves via the bristles to the Ultrasound Toothpaste medium. Millions of micro-bubbles are generated. Due to their microscopic size they are able to reach into the smallest spaces, gaps and fissures between the teeth and into the gum pockets. This also happens for wearers of implants and dental braces. The Original Ultrasound of the emmi®-dent Metallic for gentle dental cleaning and oral hygiene works without any mechanical pressure, friction or strain. Due to the powerful implosion of the micro-bubbles generated by the air oscillations of the ultrasound, all impurities such as plaque or foodstuff as well as bacteria and germs are removed efficiently. Ultrasound waves have also a positive effect on gums, gingiva and the tongue. The positive effect of Ultrasound releases also body's own healing energies, it facilitates lymphatic drainage and activates the blood flow because emmi®-dent Metallic Original Ultrasound penetrates deeply into the gum pockets. The prophylactic and healing effect of dental cleansing and oral hygiene with the Original emmi®-dent Ultrasound is generated by the effective removal of most germs and bacteria. Not only for teeth and gums, but also inside the gum tissue.

emmi®-dent Metallic cleans even where bristles and dental floss cannot. Ultrasound removes bacteria and germs deep in the gums.

Battery-charging technology

After use, the hand piece can be placed back in the charging dock for recharging the lithium battery or for a full battery charge for 12 days of toothbrush pleasure. If you use your emmi®-dent Metallic for more than 6 minutes per day, e.g. in case of several users, take a shorter interval for recharging the battery. If the battery is nearly empty, the power indicator does not light up when switching on the emmi®-dent Metallic (bottom green LED). In this case the red LED shortly flashes for 5 seconds and the device is switched off completely. After recharging the device for 24 hours, the normal functions are restored.

Battery exchange and disposal

The unit contains a special, environment-friendly, long-lasting battery. The battery is contained in the handpiece of the device. If after long-term use of the unit, the red LED flashes in spite of extensive recharging, the battery must be exchanged. Only use an emmi®-dent Metallic Lithium-Save-battery for the emmi®-dent Metallic (original exchange part with internal safety switch). If you cannot exchange the battery yourself, we will be happy to do this for you in a professional manner with low costs. If you exchange the battery yourself, please do the following: Loosen the screw at the bottom of the handpiece and remove the battery holder jointly with the bottom lid. When inserting a new battery, check the correct polarity. Insert the positive terminal in the bottom lid holder. During assembly, please ensure that the parts are correctly seated and sealed. Only then, safety class IP54 is guaranteed. You can obtain replacement batteries from Emmi Ultrasonic GmbH or in the relevant trade.

Note: The warranty, if still applicable, is only valid when the product is handled properly and the original spare parts are used.

The special battery in the device may not be disposed as normal domestic waste. For disposal, community collection points as well as other containers are available (e.g. at all sale points for similar batteries). Please see the hazard instructions mentioned below.

Hazard instruction

The emmi®-dent Metallic contains a battery with build-in safety technology. Even then, the battery may not come to kids' hands and the instructions for disposal must be followed. If in spite of all safety measures batteries are leaking, avoid contact with your skin and dispose the battery as described above. Safety measures are also required for the connection to the power supply (socket). The charging dock has a fully insulated cast housing, but never immerse the charging dock in water. The connecting cable must be in flawless condition. In case of damage to the housing, plug or connections, the whole charging dock must be discarded and replaced. Only use the emmi®-dent Metallic for its intended use as described. Children may only use the device according to the instructions and with the supervision of an adult.

Warranty

The emmi®-dent Metallic has a warranty of 24 months from the date of purchase, and a warranty of 6 months for the battery from the date of purchase. Any defects in the material and/or manufacturing errors that occur during this warranty term will be resolved without any costs. Your unit will be repaired or replaced. To the extent permitted by law, all other claims are excluded. The warranty will not apply in case of interference or dismantling by third parties without the prior approval of Emmi Ultrasonic GmbH. The warranty will also be void for products that have been damaged through improper handling or incorrect use. If under a warranty repair is granted, the guarantee term will not be extended. For claims that are not justified, for example, in case of installation or operating faults, we reserve the right to claim the costs incurred for the same. For claims under the warranty it is essential to provide the sales receipt stating the date of the purchase. For the fast processing of your claim, please send the complete unit including brush head and charging dock with sufficient shipping payment to Emmi Ultrasonic GmbH. Unfortunately, unfree shipments cannot be accepted. Herewith we point out the legal rights of the consumer as well as that these are not limited by the warranty. We thank you for your understanding.

Intended use

The products may only be used as described in the user manuals and technical descriptions, which is for dental hygiene. Any use of non-original parts or components is subject to the recommendations or approval by Emmi Ultrasonic GmbH. Any other use or which goes beyond that use is considered as improper use. Emmi Ultrasonic GmbH accepts no liability for any damages that may result from the non-intended use of the product, and we assume no warranty for the error-free and correct functioning of the components. Intended use also includes the observation of all instructions in the user manual. The symbol  on the type plate indicates the need to consult the operating instructions provided with the product. Trouble-free and safe operation of the product can only be ensured if the measures relating to proper transport, storage, assembly, installation and careful operation and maintenance are strictly observed. The tests according to BGV A3 must be performed according to the prescribed intervals. During the installation of the unit, ensure free access to the power socket on the charging dock.

Failure to use the product according to its intended use

Improper use of the product may result in hazard for life and limb as well as material damage. For this reason any use of the product other than according to its intended use is strictly prohibited. Emmi Ultrasonic GmbH accepts no liability for any damages that may result from improper use and assumes no warranty for the error-free and correct functioning of the components. For the proper operation of the unit, the following instructions must be followed:

- Only use the emmi®-dent Metallic for your dental cleaning and oral hygiene with the adequate voltage as described on the type plate (at the bottom of the charging dock)
- Do not use flammable or aggressive chemical cleansing products or acid solutions, etc.
- If damages due to transport are detected, do not connect the charging dock to the power supply.
- Only for indoor use at ambient temperature.

WITH 100% ULTRASOUND

Emmi Ultrasonic GmbH

Gerauer Strasse 34
D-64546 Moerfelden-Walldorf
Tel.: +49 (0) 6105 - 40 67 94
info@emmi-ultrasonic.de

date of purchase

stamp of the seller

subject to alterations

© Copyright EMAG AG. All rights reserved. Made in Germany/EU. The content of this manual is subject to copyright. The operating instructions may not be copied in whole or in part. Changes, abbreviations, extensions and additions as well as any publication or translation by third parties require the prior written consent of EMAG AG.

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© Copyright EMAG AG. Alle Rechte vorbehalten. Hergestellt in Deutschland/ EU. Der Inhalt dieser Bedienungsanleitung unterliegt dem Urheberrecht. Die Bedienungsanleitung darf weder ganz noch teilweise kopiert werden. Veränderungen, Kürzungen, Erweiterungen und Ergänzungen sowie jede Veröffentlichung oder Übersetzung durch Dritte bedarf der vorherigen schriftlichen Einwilligung durch die EMAG AG.